

Fully automated, high-capacity hemostasis analyzer to manage the diagnostic and monitoring requirements of clinical laboratories with mid to high workload



Yumizen

G1550




Accurate Results, efficient Workflow

Safe, easy to handle and reliable for greater efficiency

Yumizen
G1550


The Yumizen G1550 automated analyzer meets those daily expectations thanks to its unique features:

- Clotting, immunologic and chromogenic assays
- Integrated nephelometric and turbidimetric measuring channels
- Up to 8 measuring channels with 2 wavelength combination
- Continuous loading (sample racks, reagents, cuvettes)
- Only one cuvette per test
- Real-time remote access
- Comprehensive system (reagents, QC, consumables, services support)
- Safety single use cuvette
- Touch screen
- Full quality record for traceability (reagents lot, calibration, QC, rerun, users login ...)

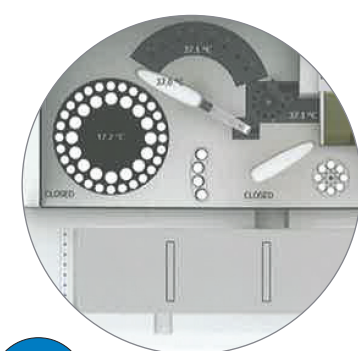

System auto-verification

- . 44 cooled controlled positions for reagents
- . Up to 4 buffer and 8 control positions
- . Liquid level sensor
- . Positive barcode identification
- . On board reagents, consumables and liquid waste management
- . Up to 1000 cuvettes for testing and 4,5L of cleaner

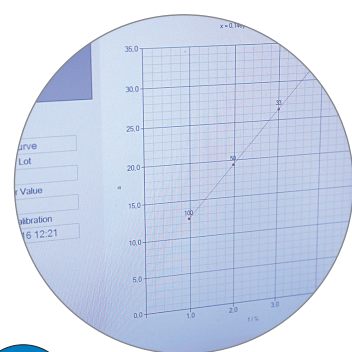


Run Calibration / QC

- . Pre calibrated tests
- . Auto-control setup, customizable according to activity (numbers of tests, new reagents ...)
- . Levey-Jennings quality control chart




Load Samples

- . 150 samples/ 15 racks
- . Cap piercing (optional) for the Closed-tube mode
- . 8 STAT sample positions

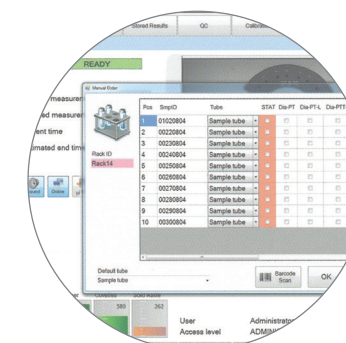

Request tests

- . Worklist
- . Random Access
- . Comprehensive test panel

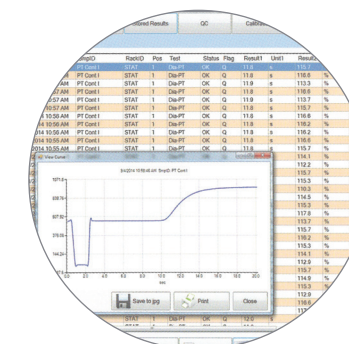

Run and validation

- . Real time curve monitoring
- . Automated rerun and reflex testing
- . Measuring status feedback

Wide range of reagents:
with pre calibrated routine tests and specialized tests (PS Ac, PC Ac, AT, DDi, Lupus Anticoagulants...)



Sample tubes rack
management for
most citrated tubes
(3.2% - 1:10)



Optional
cap-piercing*



By **HORIBA**
Medical

Hemostasis Yumizen G range



Yumizen G range

Staying fluid in Hemostasis

Dimensions & Weight:

Width : 1000 mm, Height : 600, Depth : 900 mm, Weight : 85 kg

External printer

Available

Throughput:

Yumizen G1550: 120 PT/APTT per hour

Yumizen G1500: 170 PT per hour

Sampling volume : 5-200µl

Temperature:

15°C-32°C

Power:

100-240v 50Hz-60Hz, 5A

Power supply: external power supply

Export:

6 x USB plug

Screen:

Color LCD touch screen 270x340 mm

Storage capacity:

1.000.000 results

Parameter:

Clotting tests (640 nm)

Screening tests (PT, APTT, Fibrinogen, TT)

Factors (II, V, VII, X, VIII, IX, XI, XII)

Turbidimetric test (570 nm)

D-Dimer

Chromogenic test (405 nm)

ATIII

Norms:

EN 61326-2-6 : 2006

IEC 61010-1 : 2001

IEC 61010-2-101 : 2002

EN ISO 13485 : 2003 and 2012

NF EN ISO 14971 - 2013

Directives:

98/79/EC (DIV)

2011/65/EU (RoHS)

1907/2006 (REACH)

2014/30/EU (Electro Magnetic Compatibility)

2014/35/EU (Voltage limits)

HORIBA
Medical

FRANCE +33 (0)4 67 14 15 15 - ITALY +39 / 06 51 59 22 1 - SPAIN +34 / 91- 353 30 10 - PORTUGAL +351 / 2 14 72 17 70 - UK +44 (0) 1604 542650

POLAND +48 / 22 30 05 460 - USA +1 / 949 453 0500 - BRAZIL +55 / 11 2923-5439 - THAILAND +66 / 2 861 59 95 - CHINA +86 / 21 3222 1818

INDIA +91 / 11 4646 5000 - GERMANY AXON LAB AG +49 / 7153 92260 - DISTRIBUTORS NETWORK +33 (0)4 67 14 15 16

HORIBA Medical online : <http://www.horiba.com/medical>



Operating IMS

